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QCUT 500 A

The Qcut 500 A is the powerhouse of the large cut-off machine family. Equipped with a highly dynamic axis drive and automatic travel cut/chop cut functions, the system cuts large components with exceptional precision, minimal contact, and outstanding speed.

- | Fast, powerful cuts: The highly dynamic axis drive allows rapid traverses of up to 80 mm/s and greatly reduces cutting time
- | Versatile applications: Cutting of particularly long samples is easily possible thanks to side openings on both sides
- | Multi-Position Process (MPP): Simultaneous positioning and automatic cutting of different components on the machine table

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SUPERIORITY IN DETAIL

1. Ergonomic control panel

- | Large 12.1" industrial touchscreen with intuitive QATM software
- | Quick setup, control of all axes with joystick

2. Robust machine body

- | Sturdy powder-coated aluminum construction
- | Integrable recirculating cooling unit
- | Automatic central lubrication of the axes

3. Large machine table

- | Made of stainless steel with integrated longitudinal and cross T-slots
- | Secure clamping of complex workpiece shapes using the T-slot table or Easy clamping system

4. Large cutting chamber

- | Large sliding door provides easy access and visibility
- | Side swing door for convenient access to cutting chamber
- | Side openings on the right and left for long samples
- | Break-resistant LED interior lighting

5. Optimized swivel arm

- | Highly efficient cooling system
- | Increased cutting capacity thanks to parallel guided cut-off wheel guard
- | Laser alignment aid for easy workpiece positioning

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STAINLESS STEEL MACHINE TABLE

- | Spacious machine table enables precise work even under heavy loads.
- | Integrated T-slots for the Easy clamping system for quick mounting of clamping tools
- | Efficient and reliable fixation of complex geometries using clamping devices
- | Load capacity up to 120 kg (including clamping devices)


LASER ALIGNMENT

- | Optional laser alignment projects the cutting line directly onto the workpiece, enabling precise positioning and exact cut placement.
- | Easy positioning of workpieces
- | Precise clamping without much effort


ROBUST MACHINE WITH LARGE CUTTING CHAMBER

- | Robust, powder-coated aluminum machine body features a large sliding door with scratch-resistant laminated safety glass and convenient side access door.
- | Machine body with integrable recirculating cooling unit and automatic central lubrication
- | Large cutting chamber access with safety interlock —

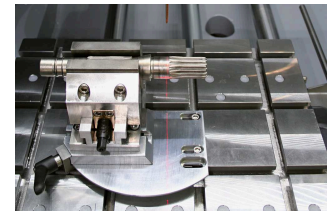PRECISION CROSS FEED

- | The table-integrated Z-axis maintains consistent cutting height while the automatic cross feed with serial cutting function delivers precise parallel cuts.
- | Simple and intuitive serial cutting
- | The multi-position process (MPP) allows more complex serial cutting to be carried out

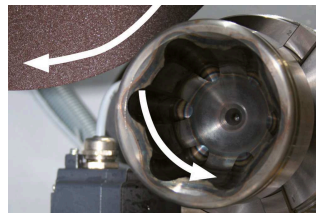ROTATION TOOL


- | The optional rotation tool features a vertical rotary feed parallel to the wheel with software-controlled angle input.
- | More uniform cutting surfaces and cutting quality
- | Cutting time is halved for rotationally symmetrical parts
- | Protects against cut-off wheel run-out

INTUITIVE AND EASY OPERATION

Large 12.1" industrial touchscreen with intuitive QATM user software and joystick.

- | Easily manage process parameters and create customer-specific maintenance tasks (200 user programs can be saved)
- | Parameters such as feed rate, travel distance, and cutting force can be changed during automatic operation

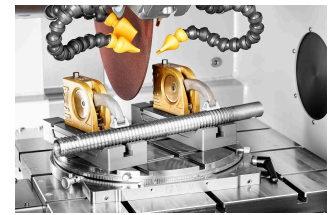PULSED CUTTING

- | For thick and difficult-to-cut materials, the adjustable automatic pulsed cutting provides gentle sectioning with optimal cooling of the cut interface.
- | Enhanced cutting results for larger components
- | Superior cooling of contact surfaces

CUTTING FORCE DEPENDENT FEED CONTROL (X- AND Y-AXIS)

The feed rate automatically decreases when cutting force reaches a preset limit, protecting the machine, sample, and clamping tools from overloading.

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
Compact, integrable cooling system for standard applications.

- | Can be integrated into the housing and extended
- | Multi-chamber filtration
- | Inlet sieve with filter bag (see consumables)

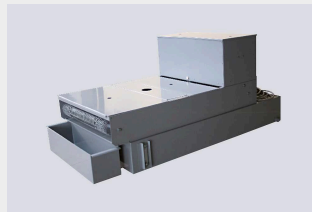BELT FILTER SYSTEM

Large cooling and filter system for continuous operation, e.g. in serial applications.

- | Tank size: 285 l
- | Filtering capacity: 100 l/min (referring emulsion)
- | Pump capacity: 100 l/min

CENTRIFUGAL SEPARATOR

Work safety: Working environment with improved air quality.

- | Removes coolant mist
- | Purifies the air and increases visibility
- | Exhaust power 452 or 500 m³/h depending on voltage

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
- | Low thermal influence on the specimen during the cutting process
- | Wet abrasive cutting achieves very fine surface finishes after cutting
- | Optimal cutting shortens the further preparation process

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- | Suitable for regular and irregularly shaped samples
- | Fast, repeat-accurate clamping of regularly shaped samples
- | Large range of different clamping tools and systems available

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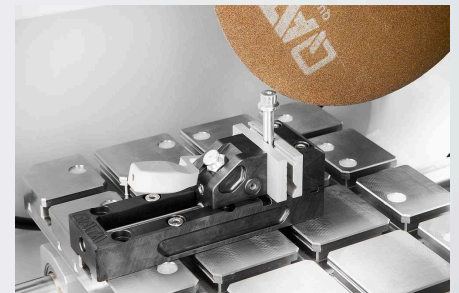QUICK CLAMPING TOOLS, CLAMPING_DEVICES

QTOOL 80, CLAMPING_DEVICES, QTOOL_SERIES

Z2231203 Qtool 80 - 100mm, Clamping_devices, Qtool_series

Z2231204 Qtool 80 - 150mm, Clamping_devices, Qtool_series

02231243 Clamping jaw 40 mm, Clamping_devices, Qtool_series

Z2231215 Clamping jaws 70 mm, Clamping_devices, Qtool_series

02231245 Clamping jaw 40 mm bevelled, Clamping_devices, Qtool_series

02231246 Clamping jaw 40 mm cross prism, Clamping_devices, Qtool_series

Z2231216 Clamping jaws 70 mm slotted, Clamping_devices, Qtool_series

Z2231217 Clamping jaws 70 mm slotted prism, Clamping_devices, Qtool_series

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Z2231211 Clamping jaws, 70 mm, Clamping_devices, Qtool_series

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Z2231212	Round bar jaws, 70 mm slotted, Clamping_devices, Qtool_series
Z2231213	Prism jaws, 70 mm slotted, Clamping_devices, Qtool_series
Z2231214	Clamping jaws, 70 mm slotted, Clamping_devices, Qtool_series

POWER CLAMPING SYSTEMS, CLAMPING_DEVICES, POWER_CLAMPING_SYSTEMS

POWER CLAMPING SYSTEM 100-540, CLAMPING_DEVICES, POWER_CLAMPING_SYSTEMS

Z2234016	Power clamping system 100-540, Clamping_devices, Power_clamping_systems
Z2234018	Prism clamping jaw, Clamping_devices, Power_clamping_systems
Z2234017	Clamping jaw flat, Clamping_devices, Power_clamping_systems
Z2234019	Angular drive, Clamping_devices, Power_clamping_systems

POWER CLAMPING SYSTEM LC 90, CLAMPING_DEVICES, POWER_CLAMPING_SYSTEMS

Z2234001	Power clamping system LC 90, Clamping_devices, Power_clamping_systems
Z2234003	Carrier jaws LC 90, Clamping_devices, Power_clamping_systems
Z2234004	V-groove jaw 100, Clamping_devices, Power_clamping_systems
Z2234005	Clamping jaw flat, Clamping_devices, Power_clamping_systems
Z2234009	Quick-change adapter, Clamping_devices, Power_clamping_systems
Z2234008	Step jaw 100, Clamping_devices, Power_clamping_systems
Z2234002	Clamping system X-Clamp 100, Clamping_devices, Power_clamping_systems
Z2234006	Socket 90°, Clamping_devices, Power_clamping_systems

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CLAMPING ACCESSORIES FOR ROTATING, SWIVELING AND SPECIAL APPLICATIONS, CLAMPING_DEVICES, CLAMPING_ACCESSORIES_SPECIALS

Z1304636	Expansion clamping device, Clamping_devices, Clamping_accessories_specials
Z1340000	クロス?バース 80-100, Clamping_devices, Clamping_accessories_specials

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ROTATIG DEVICES, CLAMPING_DEVICES, ROTATION_DEVICES

Z1870059 Rotation device type RD1, Clamping_devices, Rotation_devices

Z1891001 Rotation device type RD2.1, Clamping_devices, Rotation_devices

汎用バースとアクセサリ, CLAMPING_DEVICES, UNIVERSAL_VICES_AND_ACCESSORIES

Z1350025 Clamping rail, Clamping_devices, Universal_vices_and_accessories

01350022 Prism 120°, Clamping_devices, Universal_vices_and_accessories

01350023 Prism slim, Clamping_devices, Universal_vices_and_accessories

01350024 Prism 90°, Clamping_devices, Universal_vices_and_accessories

01350025 Prism 90° high, Clamping_devices, Universal_vices_and_accessories

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92005796 縦押さえVシュー, Clamping_devices, Universal_vices_and_accessories

92005797 横長Vシュー, Clamping_devices, Universal_vices_and_accessories

92005799 縦長ブリッジ, Clamping_devices, Universal_vices_and_accessories

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